

COLD AND HOT COMPRESS PHYSIOTHERAPISTS

一次性冰袋

Instant Cold Pack

一次性热袋

Instant Hot Pack

冷热敷理疗袋

Hot&Cold Pack




一次性冰袋 / Instant Cold Pack



一次性冰袋可以提供快速有效的冷疗，广泛应用在医院，家庭护理，急救，运动或其他户外活动。

Instant cold pack can provide quick and effective cold treatment, which is widely used in hospital, home care, first aid, sports or other outdoor activities.

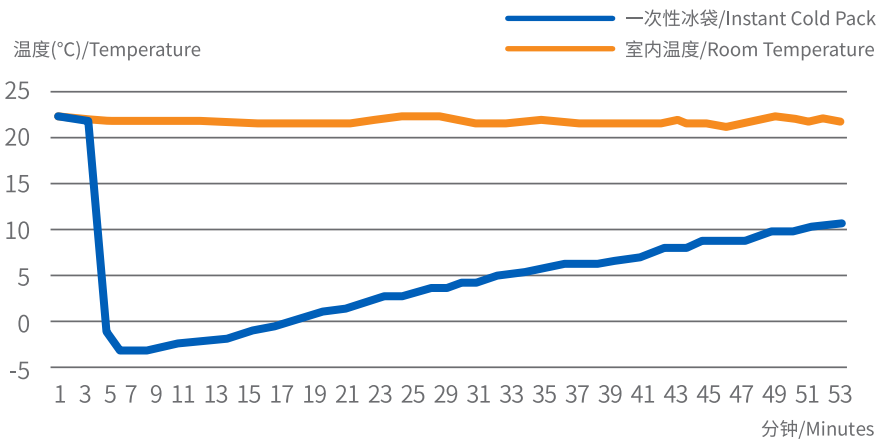

产品参数/Product parameters:



编码 Item Number	尺寸 Description	重量 Weight	Material 材质
ICE 0506	12.7*15.3CM	109G	PE&尼龙 PE&NYLON
ICE 0609	15.3*22.9CM	237G	PE&尼龙 PE&NYLON
其他/Others	可定制/Customizable		



特点/Characteristic:

用于一次性冷敷，无需制冷，即刻启动；
Single use cold application without the drippy mess of ice cubes.

仅限外部使用；无需冷藏，方便携带；
External application only. Needs no refrigeration, convenient and portable.

用于缓解由扭伤、拉伤、挫伤、轻微烧伤、牙痛和虫咬造成的疼痛和肿胀；
For relief of pain and swelling caused by sprains, strains, contusions, minor burns, toothaches and insect bites.

使用/Usage:

先用手定位到内部水袋，将里面水袋捏破；
Locate inner pouch. Squeeze firmly to rupture inner fluid pouch.

将水和内料摇晃并充分混合；
Shake briefly to mix contents.

将产品外包好护具或者毛巾后使用；
Wrap in soft cloth and apply.


在21°C~24°C左右的室温环境中使用最佳。
Best if activated about 21°C~24°C.

一次性热袋 / Instant Hot Pack



一次性热袋可以提供快速有效的热疗，广泛应用在医院，家庭护理，急救，运动或其他户外活动。

Instant hot pack can provide quick and effective hot treatment, which is widely used in hospital, home care, first aid, sports or other outdoor activities.

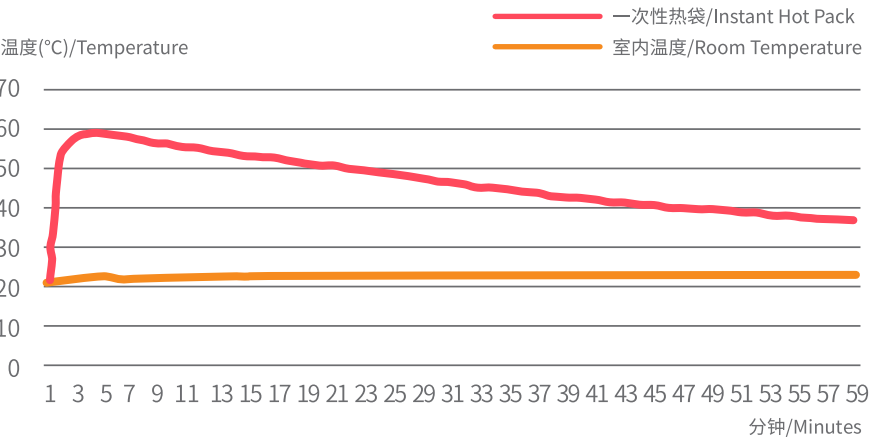

产品参数/Product parameters:



编码 Item Number	尺寸 Description	重量 Weight	材质 Material
HOT 0506	12.7*15.3CM	104G	PE&尼龙 PE&NYLON
HOT 0609	15.3*22.9CM	247G	PE&尼龙 PE&NYLON
其他/Others	可定制/Customizable		

特点/Characteristic:

无需加热，即刻启动；是用于急救的理想产品；
Start immediately without heating; it is an ideal product for first aid.

瞬间缓解撞伤、瘀伤、扭伤和运动伤；
Instant relief for knocks, bruises, sprains and sports injuries.

在寒冷/潮湿季节，可缓解由冻伤或关节炎带来的不适，改善血液循环；
Relieve frostbite in cold season, or arthritis in wet weather, Improve blood circulation.

使用/Usage:

用手取出内袋，不要用剪刀或小刀切割；
Remove the inner bag by hand, do not cut with scissor or a knife.

挤压内袋几秒钟，以激活内料；
Crush the inner bag for a few seconds to activate the contents.

把热包敷在衣服上，不要敷在皮肤上；
Apply the heat pack to clothing and not the skin.

在21°C~24°C左右的室温环境中使用最佳。
Best if activated about 21°C~24°C.

冷热敷理疗袋 / Hot & Cold Pack



冷热袋可根据身体各部位定制，冷热两用。

Hot and cold packs can be customized according to various parts of the body, and can be used for both hot and cold.

个人护理，使用方便；可重复使用；家庭必备、保健养生。

Personal care, easy to use;Can be reused for unlimited times;Family necessities, health care.

产品参数/Product parameters:



编码 Item Number	<input checked="" type="checkbox"/> HC-GEL 0510 <input type="checkbox"/> 其他/Others		
尺寸 Description	<input checked="" type="checkbox"/> 12.7*25.4CM <input type="checkbox"/> 可定制/Customizable		
重量 Weight	<input checked="" type="checkbox"/> 320G <input type="checkbox"/> 可定制/Customizable		
适用部位 Applicable parts	<input checked="" type="checkbox"/> 通用款/General Type <input type="checkbox"/> 眼部/Eye <input type="checkbox"/> 腰部/Waist <input type="checkbox"/> 关节/Joint <input type="checkbox"/> 肩颈部/Shoulder and neck		
外袋材质 material of pouch	<input checked="" type="checkbox"/> 布料/Fabric <input type="checkbox"/> PVC <input type="checkbox"/> PE&NYLON/PE尼龙		
内部填充物 Contents	<input checked="" type="checkbox"/> 凝胶/Gel  物理凝胶， 手感舒适。 The physical gel is comfortable to handle.	<input type="checkbox"/> 珠珠/Gel Beads  高分子吸水树脂，外形美观； 水量足，韧性高。 High molecular absorbent resin; Beautiful appearance; Enough water and high toughness.	<input type="checkbox"/> 粘土/Clay  天然黏土材料， 更加环保，健康。 Natural Clay, more environmental and healthier.



冷用/Cold therapy application:

提前放入冷藏或者冷冻室2小时。
Freeze within two hours.

减轻炎症引起的肿胀和疼痛，缓解由以下症状引起的不适：关节和肌肉损伤/背部和颈部疼痛、关节炎/鼻窦炎和压力性头痛、手术后疼痛/晒伤和其他烧伤。
Reduce the swelling and pain caused by inflammation and relieve the discomforts associated with:Joint & muscle injuries / Back & neck pain、Arthritis / Sinus & stress headaches、Post-surgical pain / Sunburn & other burns.

热用/Heat therapy application:

根据指定时间放入微波炉加热；如果温度不够则再加热5秒直到达到理想的温度。
Heat in microwave for specified time. Heat for additional 5 seconds as needed until desired temperature is reached.

放松紧绷的肌肉，在各种各样的疾病中帮助恢复柔韧性，例如：肌肉疼痛和酸痛、肌肉痉挛和抽筋、关节炎/月经痉挛、慢性颈部或背部僵硬。
Relax tense muscles and restore flexibility in a wide variety of ailments such as :Muscle pain & soreness、Muscle spasms & cramps、Arthritis / Menstrual cramps、Chronic neck or back stiffness.



Global INTCO Global Health

全球英科 健康全球

OUTLINE

作为冷热敷行业中首家通过FDA和CE双认证的上市企业，连续多年在全球医疗领域取得一致好评。

As the first listed enterprise in the cold and hot therapy industry to pass the certification of FDA and CE, and has achieved unanimous praise in the global medical field for many years.

全球最大的冷热敷供应商之一，医疗产品遍布全球各地。

One of the world's largest manufacturer of cold and hot pack, medical products are all over the world.

江苏英科医疗制品有限公司

Jiangsu Intco Medical Products Co., Ltd.

江苏省镇江新区烟墩山路298号
No.298 Yandunshan Road, Zhenjiang,
Jiangsu, 212132 China
电话/ Tel : 0086-511-83174088
传真/ Fax: 0086-511-83174188



英科医疗公众号



英科医疗官网